



**LIGHTHOUSE**  
WORLDWIDE SOLUTIONS  
*Technology That Counts!*



## REMOTE 5012PN

AIRBORNE PARTICLE COUNTER

4-20mA Output

The **REMOTE 5012PN** meets the stringent requirements of today's aseptic processing applications. Incorporating the latest in laser optical particle counting technology, Lighthouse designed the **REMOTE 5012PN** for continuous trouble-free operation.

The **REMOTE 5012PN** uses a built-in vacuum pump and flow control system to insure accurate and reliable sensor performance.

With a sensitivity of 0.5 micron at a controlled flow rate of 0.1 CFM, the **REMOTE 5012PN** provides continuous real-time data at a cost-effective price per point.

The **REMOTE 5012PN** integrates seamlessly into large facility monitoring and management systems and transfers up to 2 channels of simultaneous particle count data using 4-20mA output.

**Designed and built by  
Lighthouse - a name you can trust.**

Lighthouse is an ISO 9001:2000 Registered Company

[www.golighthouse.com](http://www.golighthouse.com)

### *Features:*

- ★ Channel sizes: 0.5 - 5.0 $\mu$ m for GMP Compliance
- ★ 0.1 CFM (2.83 LPM) Flow Rate
- ★ Internal Vacuum Pump with Flow Control
- ★ Up to 2 Channels of Simultaneous Data
- ★ Meets JIS Standards
- ★ NEMA Style 316L Stainless Steel Enclosure or
- ★ NEMA Style Powder Coated Steel Enclosure
- ★ LED Status Indicators
- ★ Interfaces to Facility Monitoring Systems
- ★ 4-20mA Output
- ★ 120VAC or 24VDC, 100W Models Available

### *Benefits:*

- ★ 2 Year Warranty
- ★ Cost Effective Price per Point
- ★ Easy System Integration
- ★ AC and DC Power Options
- ★ International Sales and Service Support
- ★ Low Cost of Ownership

### *Applications:*

- ★ HVAC Monitoring
- ★ Indoor Air Quality Monitoring
- ★ Food Process Monitoring
- ★ Public and Private Building Monitoring
- ★ Aseptic Process Monitoring
- ★ Isolator Monitoring
- ★ Aerospace
- ★ Chemical Manufacturing
- ★ Equipment Monitoring
- ★ Pharmaceutical Packaging and Manufacturing
- ★ Trend Analysis
- ★ Cleanroom Laundries
- ★ Cleanroom Packaging
- ★ Semiconductor
- ★ Hard Disk Drive Manufacturing



## Specifications:

<b>5012PN</b>	
<b>Size Range:</b>	0.5 - 25.0µm
<b>Channel Sizes:</b>	Standard: 0.5, 5.0µm Optional: 0.5, 10.0µm
<b>Counting Efficiency:</b>	50% @ 0.5µm, 100% for particles >0.75µm (per JIS)
<b>Flow Rate:</b>	0.1 CFM (2.83 LPM)
<b>Laser Source:</b>	Laser Diode
<b>Zero Count Level:</b>	<1 count / 5 minutes (per JIS)
<b>Concentration Limits:</b>	2,000,000 particles / ft <sup>3</sup> @ 5% coincidence loss
<b>Calibration:</b>	NIST Traceable
<b>Sample Output:</b>	Internally filtered to HEPA standards (>99.997% @ 0.3µm)
<b>Vacuum Source:</b>	Internal Pump
<b>Air Flow:</b>	Internally Monitored and Controlled
<b>Communication Modes:</b>	4-20mA Output
<b>LED Indicators:</b>	Power, Service
<b>I/O Connection:</b>	8-pin Phoenix
<b>External Pump Switch:</b>	2-pin Phoenix dry-contact relay
<b>Supporting Software:</b>	<i>Lighthouse Monitoring System; LMS XChange</i>
<b>Enclosure:</b>	NEMA Style 316L Stainless Steel or Powder Coated Steel
<b>Power:</b>	110 / 220 VAC, 50-60Hz or 24VDC 100W
<b>Dimensions:</b>	10" (w) x 12" (h) x 6.9" (d) [25.4 x 30.4 x 14.9 cm]
<b>Weight:</b>	AC Model, 20 lbs (9 kg); DC Model, 19.5 lbs (8.85 kg)

## Environmental Conditions:

**Operating:** 50°F to 104°F (10°C to 40°C), Humidity 20% to 95% non-condensing  
**Storage:** 14°F to 122°F (-10°C to 50°C), Humidity up to 98% non-condensing

## Accessories:

**Included:** Operating Manual on CD; AC or DC Power Input Connector Kit

**Optional:** Printed Operating Manual; Data Cable; Sample Tubing; 24VDC 100W Power Adapter (DC Model)


### Manufactured by:



www.golighthouse.com  
info@golighthouse.com

Tel: (800) 945 5905 (Sales and Service)  
(510) 438 0500 (Outside of USA)

### Distributed by:

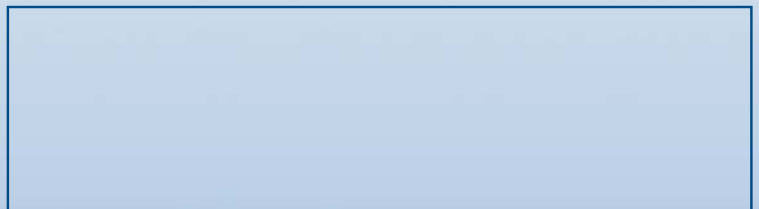

*Lighthouse Worldwide Solutions reserves the right to change specifications without notice.*

LWS PN 545401657-1 Rev 2